



A1c Go


Boronate Affinity Chromatography
POCT HbA1c Analyzer




Compact



Fully Automatic

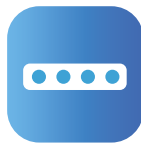

Cartridge based

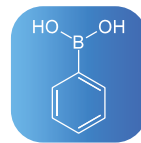


Boronate Affinity

Features

- Advanced technique: Boronate Affinity Chromatography
- Fully automatic one-step operation
- Fast and accurately answer of results in minutes
- Least interference from HbA1c Variants
- No liquid Path: no carryover & minimized maintenance
- Point-of-Care analyzer is suitable for deliver timely primary care at the physician's office

Smart Look

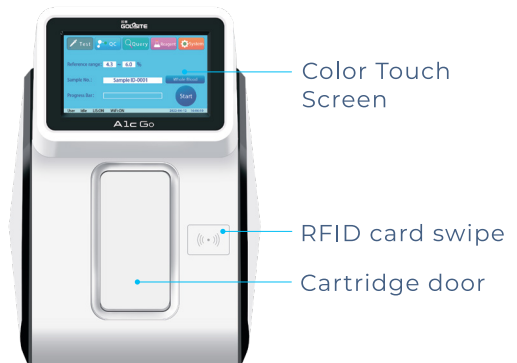

224mm(L)*160mm(W)*236mm(H)

Simplified Operation

Specifications

Screen	Color touch screen	Methodology	Boronate Affinity Chromatography
Assay	HbA1c, Glycated Hemoglobin	Pack size	25 kits/box
Sample volume	5 µL	Linear range	3.5-15.0% (NGSP)
Interval	≤5 min	Calibration	RFID card swipe
CV	≤3%	Dimensions	224mm(L)*160mm(W)*236mm(H) ≤5 Kg
Communications	LAN port, USB, Serial Port, WiFi		
consumable	Cartridge, RFID card, disposable pipette, HbA1c Control(optional)		

CONTACT

Goldsite Diagnostics Inc.
 GOLDSITE Mansion,
 No. 9 Rongtian Road, Kengzi,
 Pingshan District, Shenzhen,
 China, 518122

☎ 0086 755 2689 0807
 ✉ export@goldsite.com.cn
 🌐 en.goldsite.com.cn

@ Goldsite Diagnostics